

ENTity SD (Small Diameter) LED NasoView Nasopharyngoscope

Bringing precision endoscopy to the point of care, anywhere.



Slimline 2.7mm shaft makes the ENTity SD the preferred choice for narrow upper airway passages in both adults and children

The ENTity SD (Small Diameter) LED Nasopharyngoscope frees you to perform nasal endoscopy and laryngoscopy virtually anywhere. Our patented, integrated LED light source eliminates cumbersome external light sources while providing exceptional visibility. The result is an endoscope as portable and easy-to-use as a stethoscope. So you can take the scope to your patients, wherever they are.

New Features:

- 45X Eyepiece Magnification - *Optional*
- Slimline 2.7mm shaft



Specifications

Flexible Length	30cm (11.8")
Overall Length	53cm (20.9")
Shaft Diameter	2.7mm
Weight	11.6 oz
Articulation Range	110 degrees up / down
Depth of Field	5mm to 50mm
Field of View	85 degrees
Battery	Lithium Ion

Compact, portable, all-in-one design


- No external light source or lamp replacement
- Hybrid aluminum handle
- Fully sealed and completely immersible
- Bring the scope to the patient
- Use in the office / bedside / hospital

Better economics

- See more patients per day, increase revenue
- Reduced capital expenses and ongoing maintenance
- No lamp or light guide replacements

Superior lighting

- Patented LED light coupling
- Enables clear visibility of the true vocal chords
- Variable light intensity

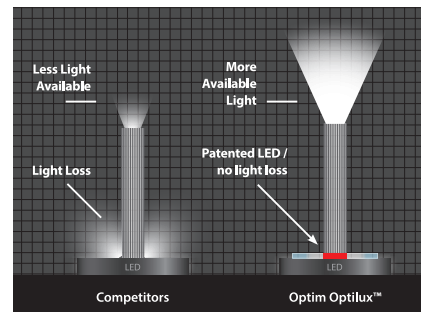

Optilux™

Optim's unique LED light source to fiber light guide coupling technology, Optilux™, in the Entity NasoView and Portable LED light source brings unmatched illumination to the end of the scope, where it counts.



64 Technology Park Road, Sturbridge, MA 01566
 Toll-free: 800.225.7486 Worldwide: 1.508.347.5100
 Fax: 508.347.2380
 Email: sales@optim-llc.com
 Website: www.optim-llc.com



Optim is ISO 9001:2008 / ISO13485:2003 Certified